

इंटरनेट

मानक

Disclosure to Promote the Right To Information

Whereas the Parliament of India has set out to provide a practical regime of right to information for citizens to secure access to information under the control of public authorities, in order to promote transparency and accountability in the working of every public authority, and whereas the attached publication of the Bureau of Indian Standards is of particular interest to the public, particularly disadvantaged communities and those engaged in the pursuit of education and knowledge, the attached public safety standard is made available to promote the timely dissemination of this information in an accurate manner to the public.

“जानने का अधिकार, जीने का अधिकार”

Mazdoor Kisan Shakti Sangathan


“The Right to Information, The Right to Live”

“पुराने को छोड़ नये के तरफ”

Jawaharlal Nehru

“Step Out From the Old to the New”

IS/ISO 10524-1 (2006): Pressure Regulators For Use With Medical Gases, Part 1: Pressure Regulators and Pressure Regulators with Flow-Metering Devices [MHD 16: Anaesthetic, Resuscitation and Allied Equipment]

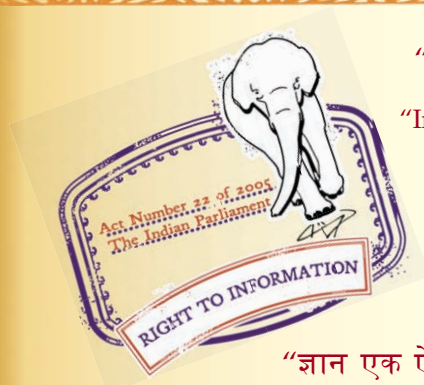

“ज्ञान से एक नये भारत का निर्माण”

Satyanarayan Gangaram Pitroda

“Invent a New India Using Knowledge”



“ज्ञान एक ऐसा खजाना है जो कभी चुराया नहीं जा सकता है”

Bhartrhari—Nitiśatakam

“Knowledge is such a treasure which cannot be stolen”

BLANK PAGE



भारतीय मानक
मैडिकल गैस के प्रेशर रेगुलेटर
भाग 1 प्रेशर रेगुलेटर एवं प्रवाह मीटरिंग यंत्रों वाले प्रेशर रेगुलेटर

Indian Standard
**PRESSURE REGULATORS FOR USE WITH
MEDICAL GASES**
**PART 1 PRESSURE REGULATORS AND PRESSURE REGULATORS
WITH FLOW-METERING DEVICES**

ICS 11.040.10

© BIS 2013
BUREAU OF INDIAN STANDARDS
MANAK BHAVAN, 9 BAHADUR SHAH ZAFAR MARG
NEW DELHI 110002

NATIONAL FOREWORD

This Indian Standard (Part 1) which is identical with ISO 10524-1 : 2006 'Pressure regulators for use with medical gases — Part 1: Pressure regulators and pressure regulators with flow-metering devices' issued by the International Organization for Standardization (ISO) was adopted by the Bureau of Indian Standards on the recommendation of the Anaesthetic, Resuscitation and Allied Equipment Sectional Committee and approval of the Medical Equipment and Hospital Planning Division Council.

This standard is published in two parts. Other part in this series is:

Part 2 Manifold and line pressure regulators

A pressure regulator is used to reduce high cylinder pressure to a lower pressure suitable for use with medical equipment or for delivery of gas directly to a patient.

These functions covers a wide range of inlet and outlet pressures and flows which require specific design characteristics. It is important that the operating characteristics of the pressure regulators are specified and tested in a defined manner.

A pressure regulator often has coupled to it, a device which controls the flow, such as a flow control valve or a fixed orifice. The flow can be indicated by a flowmeter or by a flowgauge.

It is essential that regular inspection and maintenance are undertaken to ensure that pressure regulators continue to meet the requirement of this standard.

This standard pays particular attention to:

- use of suitable materials
- safety (mechanical strength, leakage, safe relief of excess pressure and resistance to ignition)
- gas specificity
- cleanliness
- type testing
- marking
- information supplied by the manufacturer

Annex B contains rationale statements for some of the requirements of this standard. The clauses and subclauses marked with an asterisk(*) after their number have corresponding rationale, contained in Annex B, included to provide additional insight into the reasoning that led to the requirements and recommendations that have been incorporated in this standard. It is considered that knowledge of the reasons for the requirements will not only facilitate the proper application of this standard, but will expedite any subsequent revision.

The text of ISO Standard has been approved as suitable for publication as an Indian Standard without deviations. Certain conventions are, however, not identical to those used in Indian Standards. Attention is particularly drawn to the following:

- a) Wherever the words 'International Standard' appear referring to this standard, they should be read as 'Indian Standard'.
- b) Comma (,) has been used as a decimal marker while in Indian Standards, the current practice is to use a point (.) as the decimal marker.

(Continued on third cover)

Indian Standard

PRESSURE REGULATORS FOR USE WITH MEDICAL GASES

PART 1 PRESSURE REGULATORS AND PRESSURE REGULATORS WITH FLOW-METERING DEVICES

1 Scope

1.1 This part of ISO 10524 is applicable to the types of pressure regulators listed in 1.3 intended for the administration of the following medical gases in the treatment, management, diagnostic evaluation and care of patients:

- oxygen;
- nitrous oxide;
- air for breathing;
- helium;
- carbon dioxide;
- xenon;
- mixtures of the gases listed above;
- air for driving surgical tools;
- nitrogen for driving surgical tools.

1.2* These pressure regulators are intended to be fitted to cylinders with nominal filling pressures up to 25 000 kPa at 15 °C and can be provided with devices which control and measure the flow of the medical gas delivered.

1.3 The types of pressure regulators covered by this part of ISO 10524 are as follows:

- a) pressure regulators intended to be connected to cylinders by the operator;
- b) pressure regulators with integral flow-metering devices intended to be connected to cylinders by the operator;
- c) pressure regulators that are an integral part of medical equipment (e.g. anaesthetic workstations, lung ventilators, resuscitators).

2 Normative references

The following referenced documents are indispensable for the application of this document. For dated references, only the edition cited applies. For undated references, the latest edition of the referenced document (including any amendments) applies.

ISO 32:1977, *Gas cylinders for medical use — Marking for identification of content*

ISO 407:2004, *Small medical gas cylinders — Pin-index, yoke-type valve connections*

ISO 5145:2004, *Cylinder valve outlets for gases and gas mixtures — Selection and dimensioning*

ISO 5359:2000, *Low-pressure hose assemblies for use with medical gases*

ISO 9170-1:1999, *Terminal units for medical gas pipeline systems — Part 1: Terminal units for use with compressed medical gases and vacuum*

ISO 14971:2000, *Medical devices — Application of risk management to medical devices*

ISO 15001:2003, *Anaesthetic and respiratory equipment — Compatibility with oxygen*

EN 837-1:1996, *Pressure gauges — Part 1: Bourdon tube pressure gauges — Dimensions, metrology, requirements and testing*

EN 13544-2:2002, *Respiratory therapy equipment — Part 2: Tubing and connectors*

SS 01 91 02, *Colour Atlas*

3 Terms and definitions

For the purposes of this document, the following terms and definitions apply.

3.1

accuracy of flow

difference between the indicated value and the actual value of the flow expressed in percent

3.2

adjustable pressure regulator

pressure regulator that is provided with a means of operator adjustment of the outlet pressure

3.3

flow outlet

outlet intended to deliver a controlled flow of gas

3.4

flowgauge

device that measures pressure and that is calibrated in units of flow

NOTE The flowgauge does not measure flow. It indicates flow by measuring the pressure upstream of a fixed orifice.

3.5

flowmeter

device that measures and indicates the flow of a specific gas or gas mixture

3.6

gas-specific connection point

that part of the terminal unit that is the receptor for a gas-specific probe

3.7

gas-specific

having characteristics that prevent connection between different gas services

3.8

nipple

that portion of a connector that is pushed into and secured within the bore (lumen) of a hose

3.9

nominal inlet pressure

p_1

upstream pressure specified by the manufacturer for which the pressure regulator is intended to be used

NOTE For compressed gases (e.g. oxygen) p_1 is related to the cylinder filling pressure at 15 °C.

3.10

nominal outlet pressure

p_2

nominal downstream pressure

NOTE p_2 is specified by the manufacturer in the instructions for use.

3.11

orifice

restriction of known cross-section that delivers a constant flow of gas when supplied with gas at a constant upstream pressure

NOTE An orifice does not provide an indication of flow.

3.12

preset pressure regulator

pressure regulator that is not provided with a means of operator adjustment of the outlet pressure

3.13

pressure gauge

device that measures and indicates pressure

3.14

pressure outlet

outlet intended to deliver gas at a controlled pressure

3.15

pressure regulator

device that reduces the inlet pressure and maintains the set outlet pressure within specified limits

3.16

pressure-relief valve

device intended to relieve excess pressure at a preset value

3.17

single fault condition

condition in which a single means for protection against a safety hazard in equipment is defective or a single external abnormal condition is present

[IEC 60601-1]

4 Nomenclature

Examples of pressure regulators with terminology are given in Annex A.

5 General requirements

5.1 Safety

Pressure regulators shall, when transported, stored, installed, operated in normal use and maintained according to the instructions of the manufacturer, cause no safety hazard which could be foreseen using risk management procedures in accordance with ISO 14971:2000 and which is connected with their intended application, in normal condition and in single fault condition.

5.2 Alternative construction

Pressure regulators and components or parts thereof, using materials or having forms of construction different from those detailed in this clause shall be accepted if it can be demonstrated that an equivalent degree of safety is obtained.

Such evidence shall be provided by the manufacturer upon request.

NOTE Regional or national regulations can require the provision of evidence to a notified body or competent authority upon request.

5.3 Materials

5.3.1* The materials in contact with the medical gases listed in 1.1 during normal use shall be resistant to corrosion and compatible with oxygen, the other medical gases and their mixtures in the temperature range specified in 5.3.2.

NOTE 1 Corrosion resistance includes resistance against moisture and surrounding materials.

NOTE 2 Oxygen compatibility is usually defined as the ability of a material to coexist with oxygen and a moderate ignition source. The goal of using oxygen-compatible materials is to develop system designs with low probability of ignition and low consequence of ignition based on the use of materials exhibiting good compatibility and low energy release if ignited. Many materials which do not burn in air will do so in pure oxygen, particularly under pressure. Similarly, materials which can be ignited in air require lower ignition energies for ignition in oxygen. Many such materials can be ignited by friction at a valve seat or by adiabatic compression produced when oxygen at high pressure is rapidly introduced into a system initially at low pressure.

NOTE 3 Criteria for the selection of metallic and non-metallic materials are given in ISO 15001:2003.

5.3.2 The materials shall permit the pressure regulator and its components to meet the requirements of 5.4 in the temperature range of $-20\text{ }^{\circ}\text{C}$ to $+60\text{ }^{\circ}\text{C}$.

NOTE Regional or national environmental conditions might require deviation from this range of temperatures.

5.3.3 Pressure regulators shall meet the requirements of this part of ISO 10524 after being packed for transport and storage and being exposed to environmental conditions as stated by the manufacturer.

5.3.4 Springs, highly-strained components and parts liable to wear which come in contact with the medical gas shall not be plated.

NOTE Plating could come off.

5.3.5* Aluminium or aluminium alloys shall not be used for components whose surfaces come into contact with gas at cylinder pressure in normal or single-fault condition.

5.3.6 Evidence of conformity with the requirements of 5.3.1, 5.3.2, 5.3.3, 5.3.4 and 5.3.5 shall be provided by the manufacturer upon request.

NOTE Regional or national regulations might require the provision of evidence to a notified body or competent authority upon request.

5.4 Design requirements

5.4.1 Pressure gauges and flowgauges

5.4.1.1 If a Bourdon tube pressure gauge or flowgauge is used, it shall conform to EN 837-1:1996 (except for the minimum nominal size) and shall meet the requirements in 5.4.1.2, 5.4.1.3, 5.4.1.4, 5.4.1.5 and 5.4.1.6.

The requirements in 5.4.1.2, 5.4.1.3, 5.4.1.4, 5.4.1.5, 5.4.1.6 and 5.4.1.7 also apply to other types of pressure gauges and flowgauges.

5.4.1.2 If the gauge connector is threaded, it shall comply with EN 837-1:1996 or a regional or national standard.

5.4.1.3 The indicated value of a pressure gauge or flowgauge shall be legible to an operator having a visual acuity of 1 (corrected if necessary) 1 m from the gauge with an illuminance of 215 lx.

5.4.1.4 The scale of the cylinder pressure gauge shall extend to a pressure at least 33 % greater than nominal inlet pressure p_1 .

NOTE In addition to the scale ranges in EN 837-1:1996, a pressure gauge with a scale range of 0 kPa to 31 500 kPa (315 bar) can also be used.

5.4.1.5 The cylinder pressure gauge, outlet pressure gauge or flowgauge shall be class 2,5 or better in accordance with EN 837-1:1996.

5.4.1.6 The connector for a pressure gauge with a scale range greater than 4 000 kPa shall be fitted with an orifice with an area no greater than 0,1 mm².

5.4.1.7 Evidence of conformity with the requirements of 5.4.1.1 and 5.4.1.5 shall be provided by the manufacturer upon request. Compliance with the requirements of 5.4.1.2, 5.4.1.3, 5.4.1.4 and 5.4.1.6 shall be checked by visual inspection or measurement as required.

NOTE Regional or national regulations might require the provision of evidence to a notified body or competent authority upon request.

5.4.2 Connectors

5.4.2.1 Inlet connector

There shall be an inlet connector. The inlet connector for connection to cylinders shall comply with either ISO 407:1991, ISO 5145:2004 or the relevant regional or national standards. See ISO/TR 7470:1988 for information.

5.4.2.2 Outlet connector

Except for pressure regulators that are an integral part of medical equipment, the outlet connector(s) shall be in accordance with 5.4.2.2.1 and/or 5.4.2.2.2.

NOTE A pressure regulator can have multiple outlets and can have both a pressure outlet and a flow outlet.

5.4.2.2.1* Flow outlet

A flow outlet shall be fitted with a fixed nipple or a threaded connector.

Nipples, if used, shall be in accordance with EN 13544-2:2002.

Threaded connectors for oxygen or air for breathing, if used, shall be in accordance with EN 13544-2:2002 or equivalent regional or national standards or shall be manufacturer-specific connectors. Threaded connectors, if used for other gases, shall be in accordance with regional or national standards or shall be manufacturer-specific connectors.

A flow outlet shall not be fitted on a pressure regulator intended for use with air or nitrogen for driving surgical tools.

5.4.2.2.2 Pressure outlet

A pressure outlet shall be fitted with one of the following:

- a) a terminal unit or a gas-specific connection point in accordance with ISO 9170-1:1999, for the following medical gases:
 - oxygen;
 - nitrous oxide;
 - air for breathing;
 - carbon dioxide;
 - oxygen/nitrous oxide mixture 50/50 (volume fraction, %);
 - air for driving surgical tools;
 - nitrogen for driving surgical tools;
 - other gases for which terminal units in national standards exist (the connection of the terminal unit or the gas-specific connection point to the pressure regulator body need not be gas-specific).
- b) a NIST or DISS body in accordance with ISO 5359:2000, unless a national standard exists for terminal units, for the following medical gases:
 - helium;
 - xenon;
 - mixtures of oxygen and nitrous oxide [except 50/50 (volume fraction, %)];
 - mixtures of oxygen and helium;
 - mixtures of oxygen and carbon dioxide;
- c) a connector in accordance with regional or national standards.

5.4.3* Outlet pressure

Except for pressure regulators that are an integral part of medical equipment, the pressure requirements for a pressure outlet are given in 5.4.3.1.1 and 5.4.3.1.2.

The pressure requirement for a flow outlet is given in 5.4.3.2.

5.4.3.1 Pressure outlet

If a pressure regulator is fitted with a pressure outlet, the pressure regulator shall be preset.

5.4.3.1.1 Nominal outlet pressure (p_2)

Except for pressure regulators that are an integral part of medical equipment, the range of nominal outlet pressure p_2 shall be:

- a) $(400 \begin{smallmatrix} +100 \\ 0 \end{smallmatrix})$ kPa for medical gases other than air or nitrogen for driving surgical tools, or
- b) $(800 \begin{smallmatrix} +200 \\ -100 \end{smallmatrix})$ kPa for air or nitrogen for driving surgical tools.

5.4.3.1.2 Outlet pressure limits

The outlet pressure from a pressure regulator fitted with a pressure outlet (except for air or nitrogen for driving surgical tools) shall not be less than 360 kPa and not greater than 550 kPa at any flow between 0 l/min and 40 l/min for all inlet pressures between p_1 and 1 000 kPa.

The outlet pressure for a pressure regulator for air or nitrogen for driving surgical tools fitted with a pressure outlet shall not be less than 595 kPa and not greater than 1 150 kPa at any flow between 0 l/min and 350 l/min for all inlet pressures between p_1 and 2 000 kPa.

On a pressure regulator fitted with multiple pressure outlets, each pressure outlet shall be capable of meeting these requirements while all outlets are operating simultaneously.

The test for outlet pressure limits is given in 6.2.2.

5.4.3.2 Flow outlet

The pressure immediately upstream of a flow control device shall not be greater than 550 kPa for inlet pressures between p_1 and 1 000 kPa for all flow settings including zero flow.

The test for the flow outlet pressure limit is given in 6.2.3.

5.4.4 Cylinder pressure or content indicator

The pressure regulator shall be fitted with, or connected to, a cylinder pressure gauge or with an equivalent means to indicate the cylinder gas pressure or content.

NOTE In a cylinder with liquefiable gas (e.g. nitrous oxide) the pressure might not indicate the volume of gas in the cylinder.

5.4.5 Flow-metering device

If the pressure regulator is fitted with a flow outlet(s) in accordance with 5.4.2.2.1 it shall also be fitted with a flow-metering device.

A typical flow-metering device is one of the following:

- a) a flowmeter and a flow control valve supplied by a preset pressure regulator (see 5.4.12);
- b) a flowgauge and a fixed orifice supplied by an adjustable pressure regulator (see 5.4.13);
- c) multiple fixed orifices, with a means of selecting an orifice, supplied by a preset pressure regulator (see 5.4.14).

5.4.6 Pressure-adjusting device

5.4.6.1 If a pressure-adjusting device is fitted, it shall be captive such that it cannot be disengaged without the use of a tool.

Compliance shall be tested by attempting to remove the pressure-adjusting device without the use of a tool.

5.4.6.2 The pressure regulator shall be designed so that the pressure regulator valve cannot be held in the open position as a consequence of the pressure regulator spring being compressed to its solid length.

Compliance shall be verified by functional testing

5.4.6.3 Using the pressure-adjusting device it shall not be possible to set a pressure at which the pressure-relief valve opens.

Compliance shall be verified by functional testing.

5.4.7* Filtration

The pressure regulator shall be fitted on the inlet side with a filter that prevents particles greater than 100 µm from entering the pressure regulator.

NOTE 1 The selection of the material of the filter is particularly important to reduce the probability of ignition. See ISO 15001:2003 for guidance.

Evidence of conformity shall be provided by the manufacturer upon request.

NOTE 2 Regional or national regulations may require the provision of evidence to a notified body or competent authority upon request.

5.4.8* Pressure-relief valve

A pressure-relief valve shall be provided as a component part of the pressure regulator. The pressure-relief valve shall be either:

- preset; or
- not adjustable without the use of a special tool.

The leakage from the pressure-relief valve shall comply with the requirement of 5.4.9.1 up to a pressure of 550 kPa for compressed medical gases (except for air or nitrogen for driving surgical tools) and 1 150 kPa for air or nitrogen for driving surgical tools.

The pressure-relief valve shall lift automatically to relieve excess pressure and shall reset at a pressure equal to or above 550 kPa for compressed medical gases (except for air or nitrogen for driving surgical tools) and 1 150 kPa for air or nitrogen for driving surgical tools.

The discharge from the pressure-relief valve shall be equal to or greater than the maximum predicted flow through the pressure regulator valve in single fault condition at an outlet pressure of 1 000 kPa for compressed medical gases (except for air or nitrogen for driving surgical tools) and 2 000 kPa for air or nitrogen for driving surgical tools.

The pressure-relief valve shall be fitted in such a way that gas will be discharged safely.

NOTE 1 For guidance on safety, see ISO 14971:2000.

The maximum predicted flow through the pressure regulator valve in single fault condition shall be determined by the manufacturer and made available upon request.

NOTE 2 Possible single fault conditions that could result in maximum predicted flow include damage to, or loss of, the valve soft seat material.

NOTE 3 The performance characteristics of the pressure-relief valve on pressure regulators that are an integral part of medical equipment are at the discretion of the medical equipment manufacturer.

The test for the pressure-relief valve is given in 6.3.

5.4.9 Leakage

5.4.9.1 The total external leakage to atmosphere shall not exceed 0,2 ml/min (equivalent to a pressure decay of 0,020 2 kPa·l/min).

The test for total external leakage is given in 6.4.1.

5.4.9.2 The internal leakage through the pressure regulator valve shall not exceed 0,2 ml/min (equivalent to a pressure decay of 0,020 2 kPa·l/min).

The test for internal leakage is given in 6.4.2.

5.4.10 Mechanical strength

5.4.10.1 The high-pressure side of the pressure regulator shall be capable of withstanding $2,25 \times$ its nominal inlet pressure (p_1) for 5 min without rupturing.

The test is given in 6.5.1.

5.4.10.2 The low-pressure side of the pressure regulator including any integral flow control device (except for air or nitrogen for driving surgical tools) shall be capable of withstanding a pressure of 2 200 kPa for 5 min without rupturing.

The test is given in 6.5.2.1.

NOTE 2 200 kPa is $4 \times$ the maximum permissible outlet pressure of 550 kPa (see 5.4.3.1.2).

5.4.10.3 The low-pressure side of the pressure regulator for air or nitrogen for driving surgical tools, shall be capable of withstanding a pressure of 4 600 kPa for 5 min without rupturing.

The test is given in 6.5.2.1.

NOTE 4 600 kPa is $4 \times$ the maximum permissible outlet pressure of 1 150 kPa (see 5.4.3.1.2).

5.4.10.4 Components of the pressure regulator shall not be ejected if the low-pressure chamber of the pressure regulator is exposed to nominal inlet pressure (p_1) (e.g. if the regulator valve is held in the open position and the outlet connector is closed). The high-pressure gas shall either be safely retained or vented.

The test is given in 6.5.2.2.

5.4.11* Resistance to ignition

Pressure regulators for all medical gases shall not ignite or show internal scorching when subjected to oxygen pressure shocks.

The tests for resistance to ignition are given in 6.6.

5.4.12 Requirements for pressure regulators with integral flowmeters

5.4.12.1 Scales and indicators of flowmeters

Flowmeters shall be graduated in units of litres per minute (l/min) or for flows equal to or less than 1 l/min in units of millilitres per minute (ml/min).

Compliance shall be checked by visual inspection.

5.4.12.2 Legibility

The indicated value of the flowmeter shall be legible to an operator having visual acuity of 1 (corrected if necessary), 1 m from the flowmeter with an illuminance of 215 lx.

5.4.12.3 Accuracy of flow

The accuracy of the flow at any graduation of a flowmeter shall be within $\pm 10\%$ of the indicated value for flows between 10 % and 100 % of full scale or $\pm 0,5$ l/min, whichever is greater, when the flow is discharged into ambient atmosphere and corrected to reference conditions (see 6.1.3).

The accuracy of the flow at any flow graduation of a flowmeter with a maximum flow of 1 l/min or less shall be within $\pm 10\%$ of full scale.

The test for accuracy of flow is given in 6.7.

To enhance accuracy and to reduce the hazard of electrostatic discharge, means should be provided to minimize the build-up of electrostatic charges both inside and outside the flowmeter tube and its housing.

5.4.12.4 Stability of flow

The actual flow, at the maximum flow specified by the manufacturer, shall not vary by more than $\pm 20\%$ with the inlet pressure decreasing from p_1 to 1 000 kPa.

The test for stability of flow is given in 6.8.

5.4.12.5 Flow control valve

5.4.12.5.1 If a flow control valve is fitted, the flow control knob and the valve spindle shall be captive such that they cannot be disengaged without the use of a tool.

Compliance shall be tested by attempting to remove the knob and spindle without the use of a tool.

5.4.12.5.2 The flow control valve shall be designed so that the flow increases when the knob is turned anticlockwise.

Compliance shall be checked by visual inspection.

5.4.13 Requirements for pressure regulators fitted with flowgauges

5.4.13.1 Calibration

The flowgauge shall be calibrated for the identified fixed orifice and graduated in units of litres per minute (l/min).

5.4.13.2 Accuracy of flow

The accuracy of the flow at any graduation of a flowgauge shall be within $\pm 10\%$ of the indicated value for flows between 10 % and 100 % of full scale or $\pm 0,5$ l/min, whichever is greater, when the flow is discharged into ambient atmosphere and corrected to reference conditions (see 6.1.3).

The test for accuracy of flow is given in 6.7.

5.4.13.3 Stability of flow

The actual flow, at the maximum flow specified by the manufacturer, shall not vary by more than $\pm 20\%$ with the inlet pressure decreasing from p_1 to 1 000 kPa.

The test for stability of flow is given in 6.8.

5.4.14 Requirements for pressure regulators fitted with fixed orifices

5.4.14.1 Stability and accuracy of flow

The actual flow shall be within $\pm 20\%$ of each stated value or $\pm 30\%$ of each stated value for flows of 1,5 l/min or less with the inlet pressure decreasing from p_1 to 1 000 kPa.

The test for stability and accuracy of flow is given in 6.9.

5.4.14.2* Flow setting torque

If there are multiple orifices, the tangential force required at the maximum radius of the flow-selecting device to change from the "off" position and from one setting to another shall be not less than 5 N and not more than 50 N.

The test for flow setting torque is given in 6.10.

It is recommended that the flow-selecting device be designed to self-centre on a flow setting and to minimise the likelihood of selection of positions of no flow (e.g. between adjacent settings) except for the zero flow setting.

5.4.14.3 Removal of a fixed orifice

Removal of a fixed orifice shall require the use of a tool.

Compliance shall be tested by attempting to remove a fixed orifice without the use of a tool.

5.4.14.4 Legibility

The set value of the fixed orifice in use shall be legible to an operator having visual acuity of 1 (corrected if necessary), 1 m from the pressure regulator with an illuminance of 215 lx.

5.5 Constructional requirements

5.5.1* Cleanliness

Components in contact with the medical gases during normal use of pressure regulators for all gases shall meet the cleanliness requirements of ISO 15001:2003.

Evidence of conformity with this requirement shall be provided by the manufacturer upon request.

NOTE Regional or national regulations might require the provision of evidence to a notified body or competent authority upon request.

5.5.2 Lubricants

If lubricants are used, they shall be compatible with oxygen, the other medical gases listed in 1.1 and their mixtures in the temperature range specified in 5.3.2. Evidence of conformity with this requirement shall be provided by the manufacturer upon request.

NOTE 1 Attention is drawn to Annex D of ISO 15001:2003.

NOTE 2 Regional or national regulations might require the provision of evidence to a notified body or competent authority upon request.

5.5.3 Loosening torques

5.5.3.1 The torque required to remove the inlet connector from the pressure regulator body shall be $\geq 35 \text{ N}\cdot\text{m}$.

5.5.3.2 The torque required to remove the outlet connector from the pressure regulator body shall be $\geq 12 \text{ N}\cdot\text{m}$.

5.5.3.3 The torque required to remove a flow control valve (if fitted) from the pressure regulator body shall be $\geq 20 \text{ N}\cdot\text{m}$.

5.5.3.4 The torque required to remove a pressure gauge or flowgauge from the pressure regulator body shall be $\geq 12 \text{ N}\cdot\text{m}$.

5.5.3.5 The torque required to remove the flowmeter (if fitted) from the pressure regulator body shall be $\geq 20 \text{ N}\cdot\text{m}$.

The test for loosening torques is given in 6.10.

5.5.3.6 If threaded connectors are used, the requirements given in 5.5.3.1 to 5.5.3.5 apply. If other means of connection are used, an equivalent degree of safety shall be provided.

6 Test methods

6.1 General

These tests are type tests.

6.1.1 Ambient conditions

Except where otherwise stated, carry out tests at ambient conditions.

6.1.2 Test gas

In all cases carry out tests with clean, oil-free air or nitrogen with a maximum moisture content of 50 µg/g corresponding to a dew point of –48 °C at atmospheric pressure.

When a pressure regulator is tested with a gas other than that for which it is intended, the flows shall be converted using the conversion coefficients given in Table 1.

Table 1 — Conversion coefficients

Intended gas ^a	Conversion coefficient	
	Test gas: air	Test gas: nitrogen
Air	1	0,98
Oxygen	0,95	0,93
Nitrogen	1,02	1
Nitrous oxide	0,81	0,79
Carbon dioxide	0,81	0,79
Helium	2,69	2,65
Xenon	0,47	0,46

^a Flow of intended gas = Flow of test gas × conversion coefficient.

6.1.3 Reference conditions

Correct flows to 15 °C and 101,3 kPa.

6.2 Test methods for outlet pressure

6.2.1 Test equipment

Typical test equipment is shown in Figure 1.

Ensure that all equipment, including the valve controlling the flow, has a flow capacity greater than that of the pressure regulator to be tested.

6.2.2 Test method for determining outlet pressure limits for a pressure regulator fitted with a pressure outlet

Apply pressure p_1 to the inlet of the pressure regulator for medical gases (except for air or nitrogen for driving surgical tools). Set flow at 40 l/min and measure outlet pressure immediately downstream of the pressure outlet. Reduce flow to zero and measure outlet pressure immediately downstream of the pressure outlet. Repeat the test at an inlet pressure of 1 000 kPa. Verify that all outlet pressures measured are within the limits given in 5.4.3.1.2.

Apply pressure p_1 to the inlet of the pressure regulator for air or nitrogen for driving surgical tools. Set the flow to 350 l/min and measure outlet pressure immediately downstream of the pressure outlet. Reduce the flow to zero and measure outlet pressure immediately downstream of the pressure outlet. Repeat the test at an inlet pressure of 2 000 kPa. Verify that all outlet pressures measured are within the limits given in 5.4.3.1.2.

6.2.3 Test method for determining the outlet pressure limit for a pressure regulator fitted with a flow outlet

Apply a pressure of p_1 to the inlet of the pressure regulator. Set flow to maximum indicated flow and measure the pressure immediately upstream of the flow control device. Reduce the flow to zero and measure the pressure immediately upstream of the flow control device. Repeat the two pressure measurements with a pressure of 1 000 kPa applied to the inlet. Verify that all the outlet pressures are below the limit specified in 5.4.3.2.

6.3 Test method for pressure-relief valve

Apply to the low-pressure side of the pressure regulator an increasing pressure up to 550 kPa (except for air or nitrogen for driving surgical tools) or 1 150 kPa for air or nitrogen for driving surgical tools. Verify that the total leakage through the pressure-relief valve is below 0,2 ml/min (equivalent to a pressure decay of 0,020 2 kPa·l/min).

Increase the pressure until the flow reaches the maximum predicted flow (see 5.4.8). Check that the pressure complies with the requirements of 5.4.8.

Decrease the pressure to 550 kPa (except for air or nitrogen for driving surgical tools) or 1 150 kPa for air or nitrogen for driving surgical tools. Verify that at this pressure the total leakage through the pressure relief valve is below 0,2 ml/min (equivalent to a pressure decay of 0,020 2 kPa·l/min).

6.4 Test methods for leakage

6.4.1 External leakage

Measure the total external leakage of the pressure regulator at the nominal inlet pressure (p_1) and the maximum outlet pressure that can be achieved, with all outlets closed or the flow-selecting device set to zero flow. Verify that the leakage does not exceed 0,2 ml/min (equivalent to a pressure decay of 0,020 2 kPa·l/min).

6.4.2 Internal leakage

6.4.2.1 Adjustable pressure regulators

Measure the internal leakage through the pressure regulator valve at the nominal inlet pressure (p_1) with the pressure-adjusting device set to zero outlet pressure and the outlet open.

Repeat the test using an inlet pressure of 1 000 kPa for medical gases except for air or nitrogen for driving surgical tools.

Repeat the test using an inlet pressure of 2 000 kPa for air or nitrogen for driving surgical tools.

Verify that the leakage does not exceed 0,2 ml/min (equivalent to a pressure decay of 0,020 2 kPa·l/min).

6.4.2.2 Preset pressure regulators

Measure the internal leakage at the nominal inlet pressure (p_1) with the outlet plugged by monitoring the pressure.

Repeat the test using an inlet pressure of 1 000 kPa for medical gases except air or nitrogen for driving surgical tools.

Repeat the test using an inlet pressure of 2 000 kPa for air or nitrogen for driving surgical tools.

Verify that the leakage does not exceed 0,2 ml/min (equivalent to a pressure decay of 0,020 2 kPa·l/min).

6.5 Test method for mechanical strength

6.5.1 High-pressure side

For an adjustable pressure regulator, ensure that the pressure-adjusting device is in the position where the pressure regulator valve is closed.

For a preset pressure regulator, plug the outlet.

Replace the cylinder pressure gauge with a plug. Hydraulically pressurize the high-pressure side of the pressure regulator to $2,25 \times$ its nominal inlet pressure (p_1) for 5 min.

Verify that the requirements of 5.4.10.1 have been met.

6.5.2 Low-pressure side

6.5.2.1 Replace the pressure-relief valve and outlet pressure gauge, if fitted, with plugs. If necessary to hold the test pressure, either provide support for the diaphragm so that it will resist the applied pressure or replace it with a sealing blind plug. Pressurize the outlet chamber of the pressure regulator to 2 200 kPa for medical gases (except for air or nitrogen for driving surgical tools) or to 4 600 kPa for air or nitrogen for driving surgical tools for 5 min. Verify that the pressure regulator has not ruptured.

6.5.2.2 The pressure regulator valve shall be held in the open position or removed and the outlet(s) of the flow outlet and/or pressure outlet(s) blanked off. Apply a pneumatic pressure of p_1 to the inlet of the pressure regulator. Verify that no components have been ejected and the gas has been safely retained or vented.

6.6 Test method for resistance to ignition

6.6.1 General

Expose a pressure regulator to pressure shocks from industrial oxygen (minimum 99,5 % purity and hydrocarbons less or equal to 10 µg/g) through the inlet connector. The test equipment is shown in Figure 2. Before starting the test the pressure regulator shall be at room temperature.

Apply a pressure shock by increasing the pressure from atmospheric to the test pressure in a time of (20 ± 5) ms measured upstream of the pressure regulator under test (at 10 in Figure 2). Use an initial test pressure of $1,2 \times$ nominal inlet pressure (p_1) at $(60 \pm 3) ^\circ\text{C}$. During the test the inlet (test) pressure shall not decrease by more than 3 %.

Apply to the pressure regulator under test a series of 20 pressure shocks at intervals of 30 s. After each pressure shock maintain the test pressure for 10 s and then bring the pressure back to atmospheric by means of the upstream outlet valve (5 in Figure 2) and hold at atmospheric pressure for at least 3 s (see Figure 3).

After the test has been completed, dismantle the pressure regulator under test and inspect all internal parts and areas for damage (e.g. evidence of ignition or scorching).

Repeat this test on two additional pressure regulators.

Verify that the requirements of 5.4.11 have been met.

NOTE This test method is derived from ISO 2503:1998.

6.6.2 Adjustable pressure regulators

Test the pressure regulator with the regulator valve open and the outlet closed.

Apply to the pressure regulator under test a second series of 20 pressure shocks as described above with the pressure regulator valve closed.

Verify that the requirements of 5.4.11 have been met.

6.6.3 Preset pressure regulators

Test the pressure regulator in the normal delivery condition (with the pressure regulator valve open) and with the outlet closed.

6.7 Test method for accuracy of flow of pressure regulators fitted with flowmeters or flowgauges

Using the equipment shown in Figure 1, at nominal inlet pressure (p_1) set the indicated flow of the flowmeter or flowgauge under test to 10 % of full scale or the lowest graduation mark. Measure the actual flow. Repeat the test at 50 % of full scale flow and at full scale flow. Verify that the measured values are within the requirements specified in 5.4.12.3 or 5.4.13.2.

6.8 Test method for the stability of flow of pressure regulators fitted with flowmeters or flowgauges

Using the equipment shown in Figure 1 with the flow control valve (7) fully open, adjust the flow to the maximum specified by the manufacturer at nominal inlet pressure (p_1). Record the flow as indicated by the flowmeter (8) with pressures p_1 , 1 000 kPa and three or more intermediate pressures. Verify that the measured values are within the requirements specified in 5.4.12.4 or 5.4.13.3.

6.9 Test method for stability and accuracy of flow of pressure regulators fitted with fixed orifices

Use the equipment described in Figure 1 with the flow control valve (7) fully open. For each fixed orifice record the flow indicated by the flowmeter (8) with pressures p_1 , 1 000 kPa and three or more intermediate pressures. Verify that the measured values are within the requirements specified in 5.4.14.1.

6.10 Test method for flow setting and loosening torques

Measure the flow setting and loosening torques using appropriate measuring devices. Verify that the requirements of 5.4.14.2 and 5.5.3 have been met.

6.11 Test method for durability of markings and colour coding

Rub markings and colour coding by hand, without undue pressure, first for 15 s with a cloth rag soaked with distilled water, then for 15 s with a cloth rag soaked with ethanol and then for 15 s with a cloth rag soaked with isopropanol. Verify that the requirements of 7.1 have been met.

7 Marking, colour coding, packaging

7.1 Marking

7.1.1 Pressure regulators and their gas-specific components shall be durably and legibly marked with the symbol of the relevant gas in accordance with Table 2. The test for the durability of markings is given in 6.11.

NOTE In addition to the symbol, the name of the gas can be used.

Table 2 — Medical gases, marking and colour coding

Name	Symbol	Colour coding ^a
Oxygen	O ₂	White ^b
Nitrous oxide	N ₂ O	Blue ^b
Air for breathing	Air ^c	Black-white ^b
Air for driving surgical tools	Air-800	Black-white ^b
Nitrogen for driving surgical tools	N ₂ -800	Black ^b
Helium	He	Brown ^b
Carbon dioxide	CO ₂	Grey ^b
Xenon	Xe	Light brown ^d
Mixtures of the above gases	e	e
^a See Annex C for national deviations for colour coding for medical gases. ^b In accordance with ISO 32:1977. ^c National languages can be used for air. ^d An example of light brown is NCS 3030-Y30 R in accordance with SS 01 91 02. ^e According to the components.		

7.1.2 In addition to the requirement of 7.1.1, the pressure regulator shall be marked with the following:

- a) the name and/or the trademark of the manufacturer or distributor (some regional regulatory authorities do not accept that identification of the distributor replaces the identification of the manufacturer);
- b) the model or type designation;
- c) means to ensure traceability such as type, batch or serial number or year of manufacture;
- d) the value of nominal inlet pressure p_1 .

For pressure regulators that are an integral part of medical equipment, the marking listed above may be on the medical equipment.

7.1.3 If a fixed orifice is designed to be removed by use of a tool, the body of the fixed orifice shall be marked with the corresponding flow in units of l/min.

7.1.4 Pressure gauges and flowgauges shall be marked with the following:

- a) means of identification (e.g. the name and/or the trademark of the manufacturer and/or distributor);
- b) the words "USE NO OIL" or the symbol (application of ISO 7000-0248) shown in Figure 4;
- c) the unit of pressure (for pressure gauges);
- d) the unit of flow (for flowgauges);
- e) the identity of the fixed orifice for which the flowgauge is calibrated.

7.1.5 Pressure-adjusting devices and flowmeter control valves (if fitted) shall be clearly and durably marked with the direction for increasing pressure or flow.

7.1.6 Compliance with 7.1.1, 7.1.2, 7.1.3, 7.1.4 and 7.1.5 shall be checked by visual inspection.

7.2 Colour coding

7.2.1 If colour coding is used, it shall be in accordance with Table 2 or regional or national standards.

NOTE Annex C shows national and regional deviations in colour coding and nomenclature for medical gases.

7.2.2 Colour coding shall be durable. The test for the durability of colour coding is given in 6.11.

7.3 Packaging

7.3.1 Pressure regulators and spare parts shall be sealed to protect against contamination and packaged to prevent damage during storage and transportation.

7.3.2 Packages shall provide a means of identification of the contents.

7.3.3 The package shall be marked with the transport and storage conditions specified by the manufacturer.

8 Information to be supplied by the manufacturer

8.1 In order to supply the necessary information for safe use, the manufacturer of the pressure regulator shall provide to his customer(s) the following information:

- a technical description;
- instructions for operation and maintenance including the recommended frequency of maintenance activities;
- a list of recommended spare parts;
- instructions for use including cleaning and disinfection procedures;
- the address of the manufacturer.


8.2 For a pressure regulator fitted with a pressure outlet, the technical description provided shall include values of nominal inlet pressure (p_1) and nominal outlet pressure (p_2).

8.3 For a pressure regulator fitted with a flow outlet, the technical description provided shall include values of nominal inlet pressure (p_1) and range of flow settings.

8.4 Instructions for operating the pressure regulator shall give detailed information needed for the safe performance including:

- functions of the controls;
- the sequence of operations and connection and disconnection of detachable parts and accessories;
- the danger of fire or explosion arising from the use of oils and greases (including hand creams, etc.) with air, oxygen and oxygen mixtures;
- the need to open and shut the cylinder valve slowly;
- a warning not to use a flow outlet for driving any medical equipment;
- a warning that, if multiple fixed orifices are fitted, no flow may be delivered if the flow-selecting device is set between adjacent settings.

Particular attention shall be given to the lubricants, if used, in the maintenance of the pressure regulator and the danger arising from changing the setting of the pressure-relief valve.

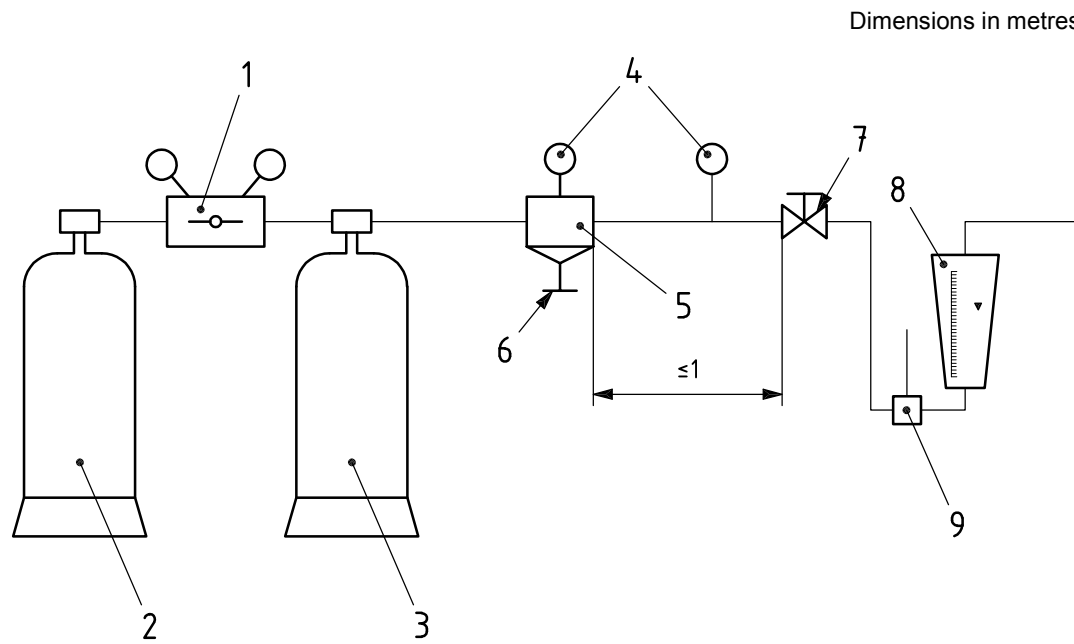


Key

- 1 auxiliary pressure regulator
- 2 gas supply
- 3 buffer cylinder
- 4 calibrated gauges
- 5 pressure regulator under test
- 6 pressure-adjusting device
- 7 flow control valve
- 8 flowmeter
- 9 thermometer

Figure 1 — Equipment for pressure and flow tests


Dimensions in millimetres

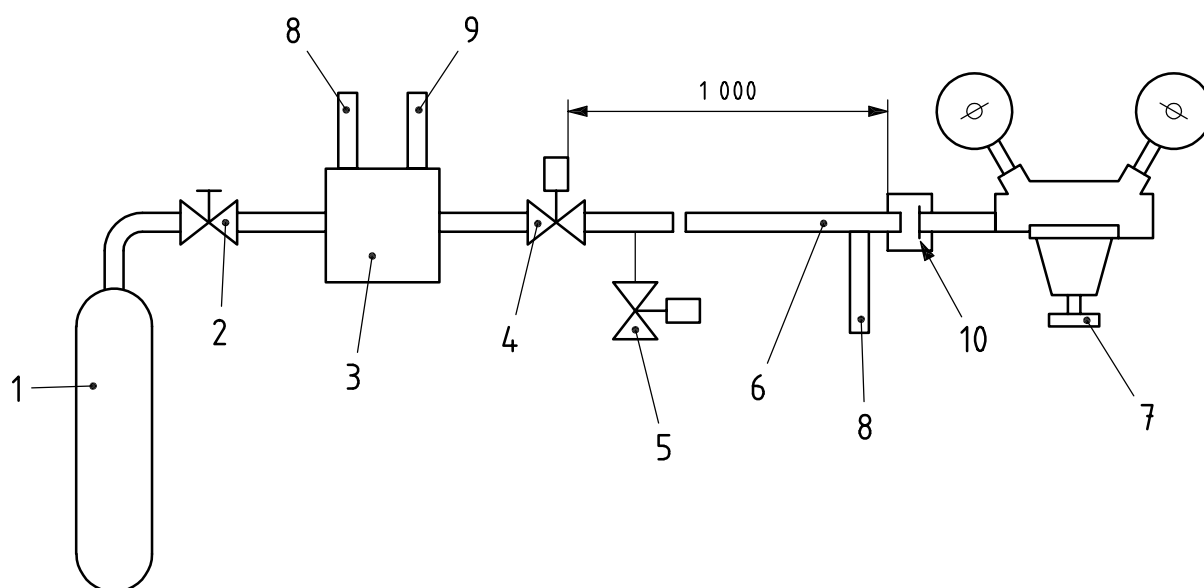

Key

- 1 oxygen supply
- 2 inlet valve
- 3 high-pressure vessel with device for pre-heating oxygen to $(60 \pm 3) ^\circ\text{C}$
- 4 quick opening valve
- 5 outlet valve
- 6 connection tube with internal diameter of 5 mm
- 7 pressure regulator under test
- 8 pressure transducer
- 9 thermometer
- 10 measuring point

Figure 2 — Test bench for test for resistance to ignition

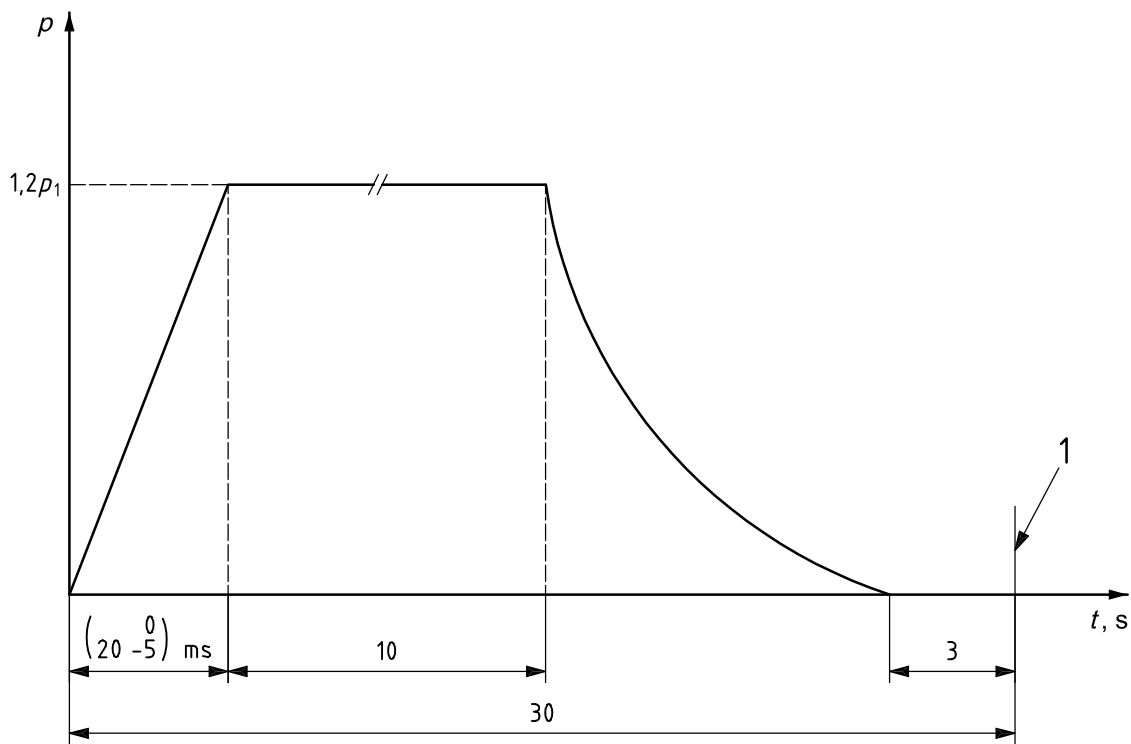

Key
 p pressure
 t time, s
1 next pressure shock

Figure 3 — Test interval


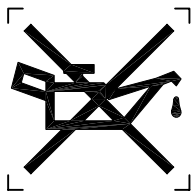


Figure 4 — Symbol for use no oil (application of ISO 7000-0248)

Annex A (informative)


Typical examples of pressure regulators and pressure regulators with flow-metering devices

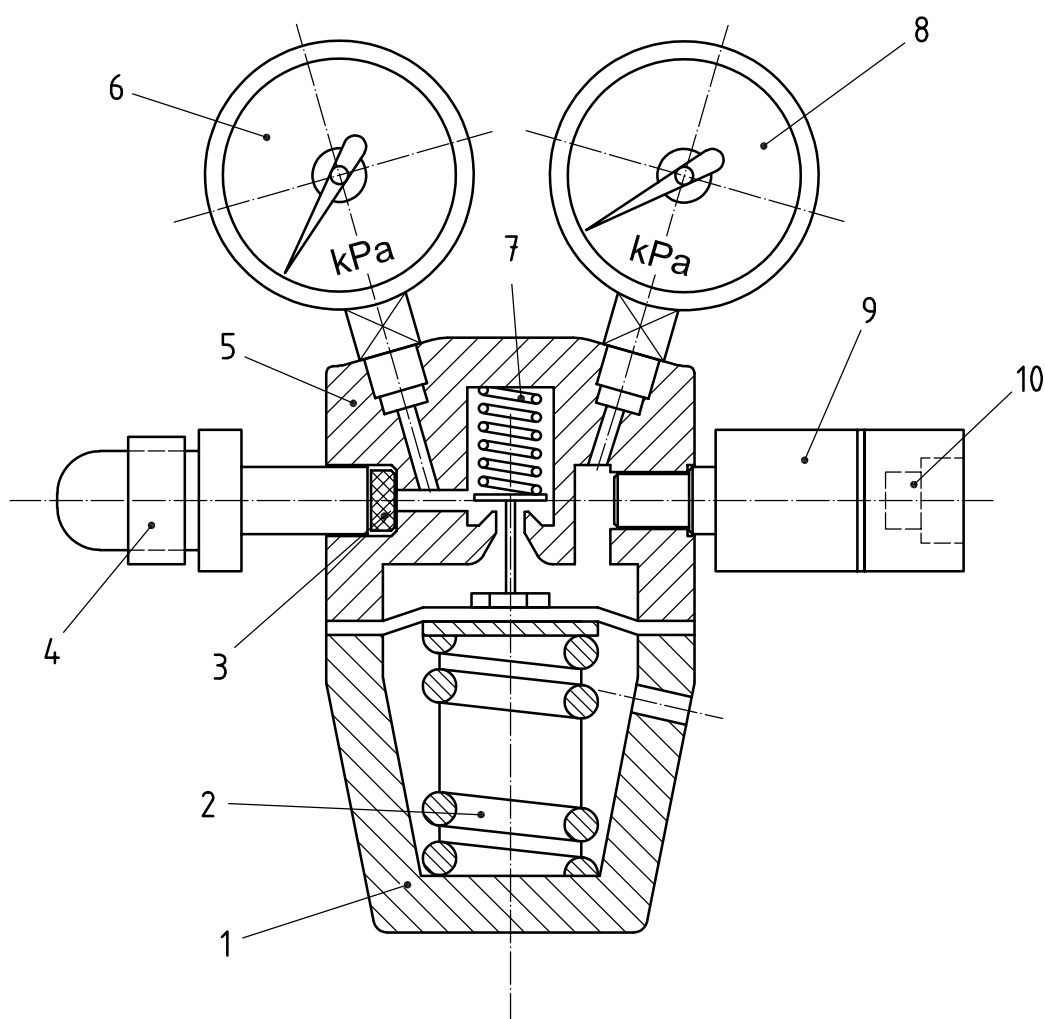

Key


- 1 cover
- 2 pressure regulator spring
- 3 inlet filter
- 4 inlet connector
- 5 body
- 6 high-pressure gauge
- 7 pressure regulator valve spring
- 8 low-pressure gauge
- 9 terminal unit
- 10 gas-specific connection point

Figure A.1 — Preset pressure regulator with terminal unit

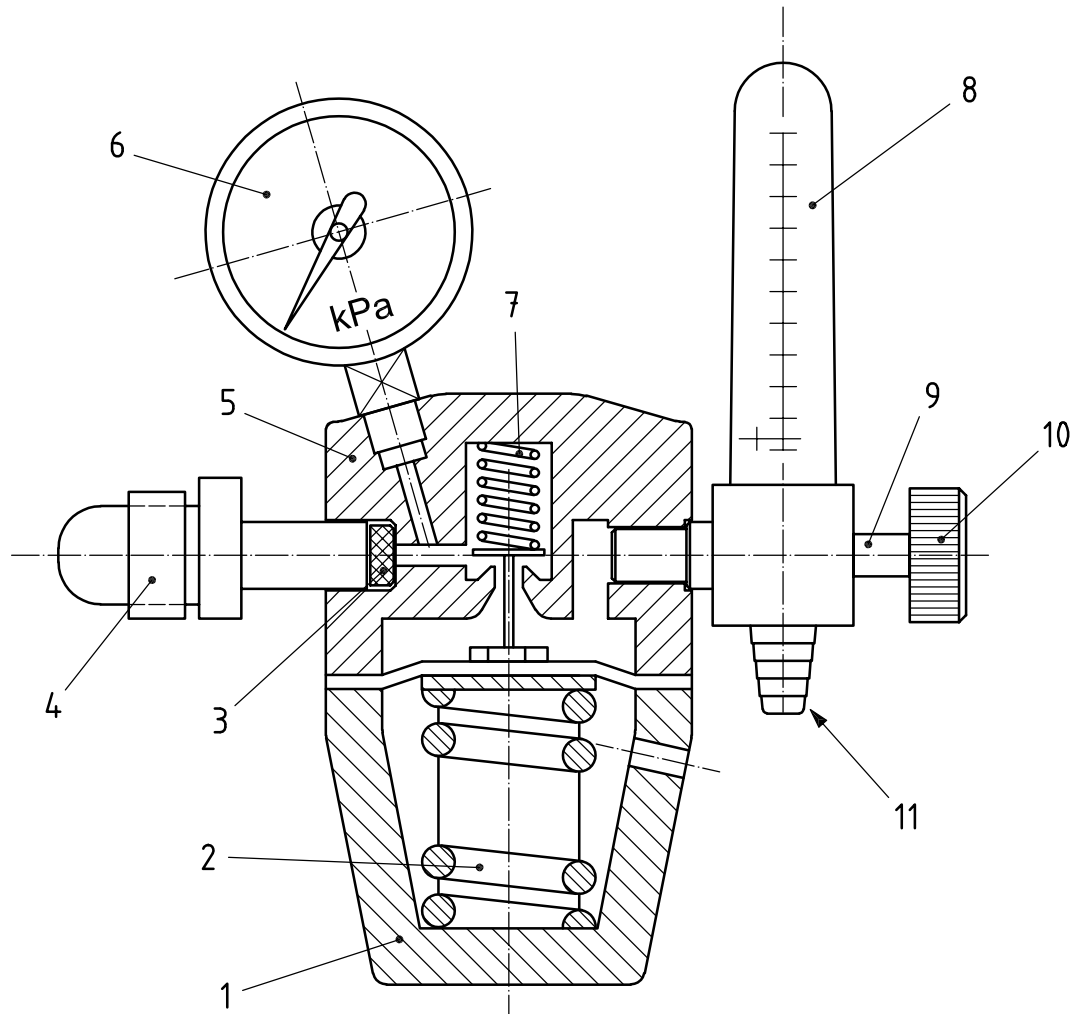**Key**

- 1 cover
- 2 pressure regulator spring
- 3 inlet filter
- 4 inlet connector
- 5 body
- 6 high-pressure gauge
- 7 pressure regulator valve spring
- 8 flowmeter
- 9 flow control valve spindle
- 10 flow control valve knob
- 11 nipple


Figure A.2 — Preset pressure regulator with flowmeter

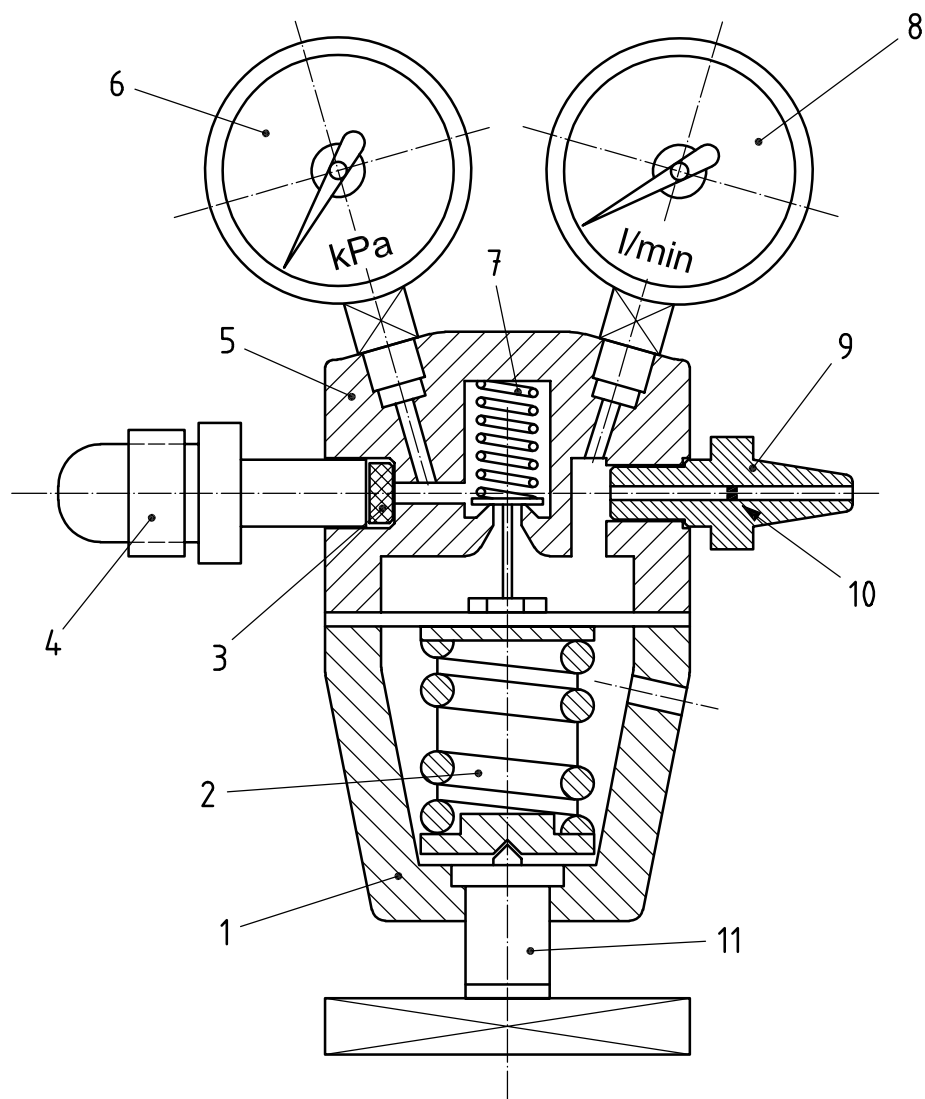

Key

- 1 cover
- 2 pressure regulator spring
- 3 inlet filter
- 4 inlet connector
- 5 body
- 6 high-pressure gauge
- 7 pressure regulator valve spring
- 8 flowgauge
- 9 body of fixed orifice
- 10 fixed orifice
- 11 pressure-adjusting device

Figure A.3 — Adjustable pressure regulator with flowgauge

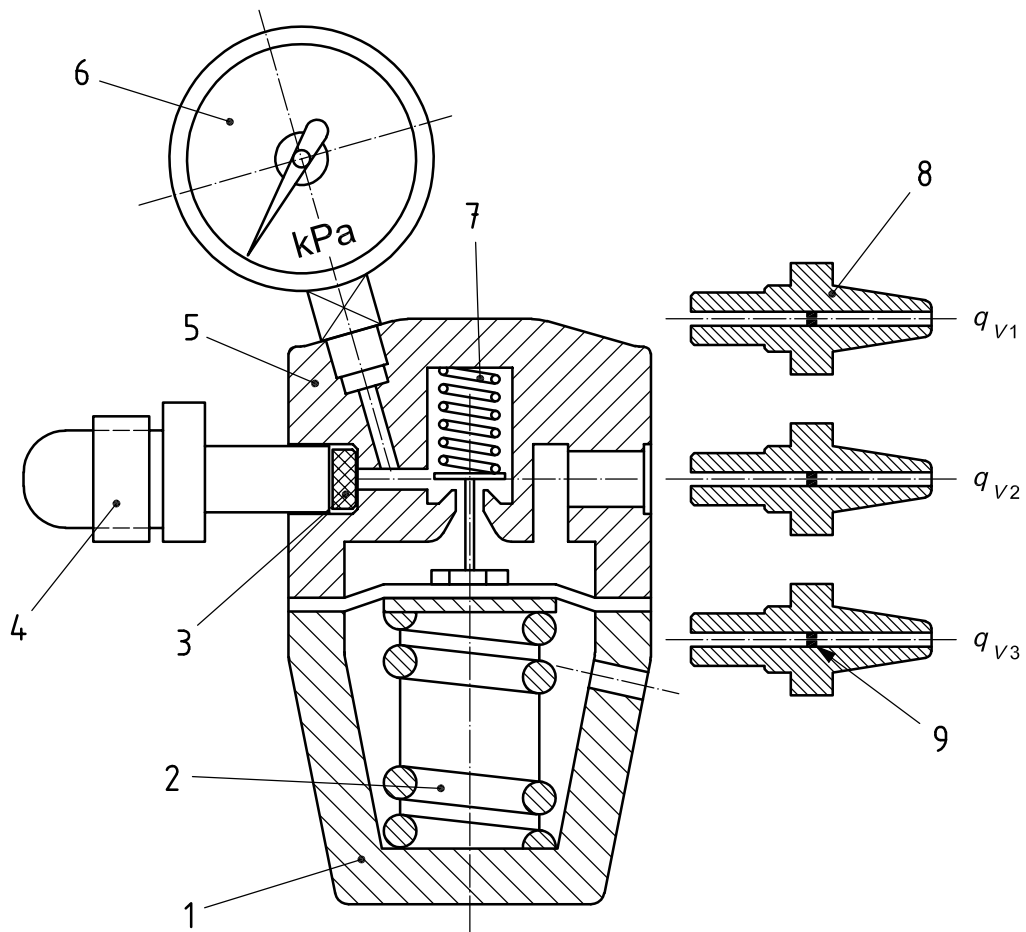**Key**

- 1 cover
- 2 pressure regulator spring
- 3 inlet filter
- 4 inlet connector
- 5 body
- 6 high-pressure gauge
- 7 pressure regulator valve spring
- 8 body of fixed orifice
- 9 fixed orifice

q_{V1} , q_{V2} , q_{V3} volume flowrates of fixed orifices, l/min.

Figure A.4 — Preset pressure regulator with multiple fixed orifices

Annex B **(informative)**

Rationale

B.1 This annex provides a rationale for some requirements of this part of ISO 10524 and is intended for those who are familiar with the subject of this part of ISO 10524 but who have not participated in its development. An understanding of the rationale underlying these requirements is considered to be essential for their proper application. Furthermore, as clinical practice and technology change, it is believed that a rationale will facilitate any revision of this part of ISO 10524 necessitated by those developments.

The following correspond to subclauses marked with * in this part of ISO 10524. The numbering is, therefore, not consecutive.

B.1.2 Cylinders used to supply medical gases are currently filled to nominal filling pressures of up to 25 000 kPa. Cylinders exist that can be filled to higher pressures (currently up to 30 000 kPa) and these are already in use for certain applications. Although these higher pressure cylinders have been used in non-medical applications, there is only limited knowledge of the requirements for their safe use. Therefore the scope of this part of ISO 10524 has been restricted to cylinders filled at pressures up to 25 000 kPa. Once experience has been gained, it is anticipated that this part of ISO 10524 will be amended to include pressure regulators for use with cylinders with nominal filling pressures of up to 30 000 kPa.

B.5.3.1 Pressure regulators for different gases are often made with interchangeable components or subassemblies. The requirement for compatibility with oxygen should therefore be applied to pressure regulators for all gases.

B.5.3.5 The bodies and parts of the high-pressure side of most pressure regulators are made of brass or aluminium. Aluminium and its alloys are more likely to ignite in an oxidizing environment than brass. In ignition tests, aluminium can burn vigorously even at low pressures, while brass burns only at pressures many times higher than cylinder filling pressures. Although there are some reported instances of ignition in brass pressure regulators, these pressure regulators have a long history of safe use and are believed to be safer than aluminium pressure regulators. Therefore components on the high-pressure side of a pressure regulator are required by this part of ISO 10524 to be composed of a material other than aluminium, e.g. brass.

Some national regulations or organizations prohibit or recommend against the use of aluminium in pressure regulators.

Pressure regulators for different gases are often made with interchangeable components or subassemblies. This requirement should therefore be applied to pressure regulators for all gases.

B.5.4.2.2.1 A flow outlet is typically used to supply a medical gas for inhalation by a patient. The flow and pressure delivered at such an outlet is not intended to drive medical equipment, which may present sufficient back pressure to make the flow outlet non-functional. Therefore a flow outlet is required to have different dimensions from a pressure outlet which is intended to drive medical equipment.

B.5.4.3 The outlet pressure has been linked to the type of outlet connector for the following reasons:

- a) When a pressure regulator is fitted with a pressure outlet, the pressure outlet should have essentially the same performance as a medical gas pipeline terminal unit. The pressure at the terminal unit is given in ISO 7396-1 which specifies the following nominal ranges:
 - 1) 400 kPa to 500 kPa with an allowable deviation of $\pm 10\%$ between conditions of zero flow and maximum flow for medical gases other than air or nitrogen for driving surgical tools;
 - 2) 700 kPa to 1 000 kPa with an allowable deviation of $\pm 15\%$ between conditions of zero flow and maximum flow for air or nitrogen for driving surgical tools.

- b) A flow outlet is not intended to supply gas to medical equipment such as a ventilator or an anaesthetic workstation. Such equipment needs to be connected to a pressure outlet.
- c) Pressure regulators fitted with NIST or DISS connectors are intended for the supply of certain medical gases which are normally not piped but which can be used for therapy or measurement.
- d) NIST or DISS connectors for those medical gases which are normally supplied by medical gas pipeline systems are not permitted by this part of ISO 10524 so that only one system for gas-specific connectors is used for any one medical gas.

B.5.4.7 Filters may be subjected to particle impact and are therefore a potential site of ignition. It is particularly important that the material of the filter itself be non-flammable under the conditions of use to reduce the likelihood of propagation of a fire. Some sintered or woven wire mesh filter materials exhibit a tendency to support combustion in laboratory tests and therefore should not be used at cylinder pressure. Therefore the selection of the filter materials should be based upon the available technical data and experimental tests.

B.5.4.8 When a pressure regulator is fitted with a pressure outlet, the pressure outlet should have essentially the same performance as a medical gas pipeline terminal unit. ISO 7396-1 specifies the following values in single fault condition:

- 1 000 kPa for compressed medical gases other than air or nitrogen for driving surgical tools (medical equipment such as ventilators and anaesthetic workstations are required to function with pressure variations up to the maximum permitted pressure);
- 2 000 kPa for air or nitrogen for driving surgical tools.

In order to avoid the application of excess pressure to downstream components, the maximum predicted flow through the pressure regulator valve in single fault condition must be known to determine the performance of the pressure relief valve.

B.5.4.11 Pressure regulators for different gases are often made with interchangeable components or subassemblies. The requirement for resistance to ignition should therefore be applied to pressure regulators for all gases.

B.5.4.14.2 A potentially dangerous situation could arise if the flow-selecting device can be unintentionally set to a position where no flow occurs. Therefore the design of the flow-selecting device should minimise the possibility of this happening. A warning of this possible hazard is required in the instructions for use.

B.5.5.1 Pressure regulators for different gases are often made with interchangeable components or subassemblies. The requirement for cleanliness should therefore be applied to pressure regulators for all gases.

Annex C (informative)

Reported regional and national deviations of colour coding and nomenclature for medical gases

Table 2 contains requirements for colour coding of medical gases in accordance with ISO 32:1977. Although many countries/markets comply with ISO 32:1977, some countries/markets have colour coding requirements that differ from those specified in ISO 32:1977. Often these alternative colour codes are mandated by standards in force within the respective countries/markets.

Table C.1 — European Union

Medical gas	Colour coding
Oxygen	White
Nitrous oxide	Blue
Medicinal air	Black and white
Nitrogen	Black
Carbon dioxide	Grey
Helium	Brown
Mixtures of gases	Combination of colours from individual gases, e.g. white/blue
NOTE See EN 1089-3 ^[10] .	

Table C.2 — United States of America

Medical gas	Colour coding
Oxygen	Green
Nitrous oxide	Blue
Medical air	Yellow
Nitrogen	Black
Carbon dioxide	Grey
Helium	Brown
Mixtures of gases	Combination of colours from individual gases, e.g. green/blue
NOTE See CGA C-9:1988 ^[5] .	

Table C.3 — Australia and New Zealand

Medical gas	Colour coding
Oxygen	White
Nitrous oxide	Ultramarine
Medical breathing air	Black and white
Surgical tool gas	Aqua
Nitrous oxide/oxygen 50/50	Ultramarine and white
Carbon dioxide	Green grey
Carbon dioxide in oxygen – nominal 5 %	White and green grey
Spare medical gas	Sand
NOTE See AS 2896-1998 ^[1] and AS 4484-2004 ^[2] .	

Table C.4 — Canada

Medical gas	Colour coding
Oxygen	White
Nitrous oxide	Blue
Medical breathing air	Black and white
Nitrogen	Black
Carbon dioxide	Grey
Helium	Brown
Mixtures of gases	Combination of colours from individual gases
NOTE See CAN/CGSB 24.2-M86 ^[4] .	

Table C.5 — Japan

Medical gas	Colour coding
Oxygen	Green
Nitrous oxide	Blue
Air for breathing	Yellow
Nitrogen	Grey
Carbon dioxide	Orange
Air for driving surgical tools	Brown
NOTE See JIS T 7101:1997 ^[15] .	

Bibliography

- [1] AS 2896-1998, *Medical gas systems — Installation and testing of non-flammable medical gas pipeline systems*
- [2] AS 4484-2004, *Gas cylinders for industrial, scientific, medical and refrigerant use — Labelling and colour coding*
- [3] ASTM G175:2003, *Standard test method for evaluating the ignition sensitivity and fault tolerance of oxygen regulators used for medical and emergency applications*
- [4] CAN/CGSB 24.2-M86, *Identification of medical gas containers, pipelines and valves*
- [5] CGA C-9:1988, *Standard color marking of compressed gas containers intended for medical use*
- [6] EN 737-1, *Medical gas pipeline systems — Part 1: Terminal units for compressed medical gases and vacuum*
- [7] EN 737-3, *Medical gas pipeline systems — Part 3: Pipelines for compressed medical gases and vacuum*
- [8] EN 738-1, *Pressure regulators for use with medical gases — Part 1: Pressure regulators and pressure regulators with flow-metering devices*
- [9] EN 739, *Low-pressure hose assemblies for use with medical gases*
- [10] EN 1089-3, *Transportable gas cylinders — Gas cylinder identification — Part 3: Colour coding*
- [11] IEC 60601-1:1998, *Medical electrical equipment — Part 1: General requirements for safety*
- [12] ISO 2503:1998, *Gas welding equipment — Pressure regulators for gas cylinders used in welding, cutting and allied processes up to 300 bar*
- [13] ISO 4135, *Anaesthetic and respiratory equipment — Vocabulary*
- [14] ISO 7396-1:2002, *Medical gas pipeline systems — Part 1: Pipelines for compressed medical gases and vacuum*
- [15] JIS T 7101:1997, *Medical gas pipeline systems*
- [16] ISO/TR 7470:1988, *Valve outlets for gas cylinders — List of provisions which are either standardized or in use*
- [17] ISO 7000, *Graphical symbols for use on equipment — Index and synopsis*

(Continued from second cover)

In this adopted standard, reference appears to certain International Standards for which Indian Standards also exist. The corresponding Indian Standards which are to be substituted in their respective places are listed below along with their degree of equivalence for the editions indicated:

<i>International Standard</i>	<i>Corresponding Indian Standard</i>	<i>Degree of Equivalence</i>
ISO 9170-1 : 1999 ¹⁾ Terminal units for medical gas pipeline systems — Part 1: Terminal units for use with compressed medical gases and vacuum	IS/ISO 9170-1 : 2008 Terminal units for medical gas pipeline systems: Part 1 Terminal units for use with compressed medical gases and vacuum	Technically Equivalent
ISO 15001 : 2003 Anaesthetic and respiratory equipment — Compatibility with oxygen	IS/ISO 15001 : 2003 Anaesthetic and respiratory equipment — Compatibility with oxygen	Identical

The technical committee has reviewed the provisions of the following International Standards/Other Publications referred in this adopted standard and has decided that they are acceptable for use in conjunction with this standard:

<i>International Standard/ Other Publication</i>	<i>Title</i>
ISO 32 : 1977	Gas cylinders for medical use — Marking for identification of content
ISO 407 : 2004	Small medical gas cylinders — Pin-index, yoke-type valve connections
ISO 5145 : 2004	Cylinder valve outlets for gases and gas mixtures — Selection and dimensioning
ISO 5359 : 2000	Low-pressure hose assemblies for use with medical gases
ISO 14971 : 2000	Medical devices — Application of risk management to medical devices
EN 837-1:1996	Pressure gauges — Part 1: Bourdon tube pressure gauges — Dimensions, metrology, requirements and testing
EN 13544-2 : 2002	Respiratory therapy equipment — Part 2: Tubing and connectors
SS 01 91 02	Colour Atlas

For the purpose of deciding whether a particular requirement of this standard is complied with, the final value, observed or calculated, expressing the result of a test, shall be rounded off in accordance with IS 2 : 1960 'Rules for rounding off numerical values (*revised*)'. The number of significant places retained in the rounded off value should be the same as that of the specified value in this standard.

¹⁾ Since revised in 2008.

Bureau of Indian Standards

BIS is a statutory institution established under the *Bureau of Indian Standards Act, 1986* to promote harmonious development of the activities of standardization, marking and quality certification of goods and attending to connected matters in the country.

Copyright

BIS has the copyright of all its publications. No part of these publications may be reproduced in any form without the prior permission in writing of BIS. This does not preclude the free use, in course of implementing the standard, of necessary details, such as symbols and sizes, type or grade designations. Enquiries relating to copyright be addressed to the Director (Publications), BIS.

Review of Indian Standards

Amendments are issued to standards as the need arises on the basis of comments. Standards are also reviewed periodically; a standard along with amendments is reaffirmed when such review indicates that no changes are needed; if the review indicates that changes are needed, it is taken up for revision. Users of Indian Standards should ascertain that they are in possession of the latest amendments or edition by referring to the latest issue of 'BIS Catalogue' and 'Standards: Monthly Additions'.

This Indian Standard has been developed from Doc No.: MHD 11 (0172).

Amendments Issued Since Publication

Amendment No.	Date of Issue	Text Affected

BUREAU OF INDIAN STANDARDS

Headquarters:

Manak Bhavan, 9 Bahadur Shah Zafar Marg, New Delhi 110002

Telephones: 2323 0131, 2323 3375, 2323 9402

Website: www.bis.org.in

Regional Offices:

Telephones

Central	: Manak Bhavan, 9 Bahadur Shah Zafar Marg NEW DELHI 110002	{ 2323 7617 2323 3841
Eastern	: 1/14, C.I.T. Scheme VII M, V.I.P. Road, Kankurgachi KOLKATA 700054	{ 2337 8499, 2337 8561 2337 8626, 2337 9120
Northern	: SCO 335-336, Sector 34-A, CHANDIGARH 160022	{ 260 3843 260 9285
Southern	: C.I.T. Campus, IV Cross Road, CHENNAI 600113	{ 2254 1216, 2254 1442 2254 2519, 2254 2315
Western	: Manakalaya, E9 MIDC, Marol, Andheri (East) MUMBAI 400093	{ 2832 9295, 2832 7858 2832 7891, 2832 7892

Branches: AHMEDABAD. BANGALORE. BHOPAL. BHUBANESHWAR. COIMBATORE. DEHRADUN. FARIDABAD. GHAZIABAD. GUWAHATI. HYDERABAD. JAIPUR. KANPUR. LUCKNOW. NAGPUR. PARWANOO. PATNA. PUNE. RAJKOT. THIRUVANATHAPURAM. VISAKHAPATNAM.